


**BRUNO**

# *Electra-Ride™ II*

*Straight Rail Stairlifts*

**B R E A K I N G   T H E   B A R R I E R S**

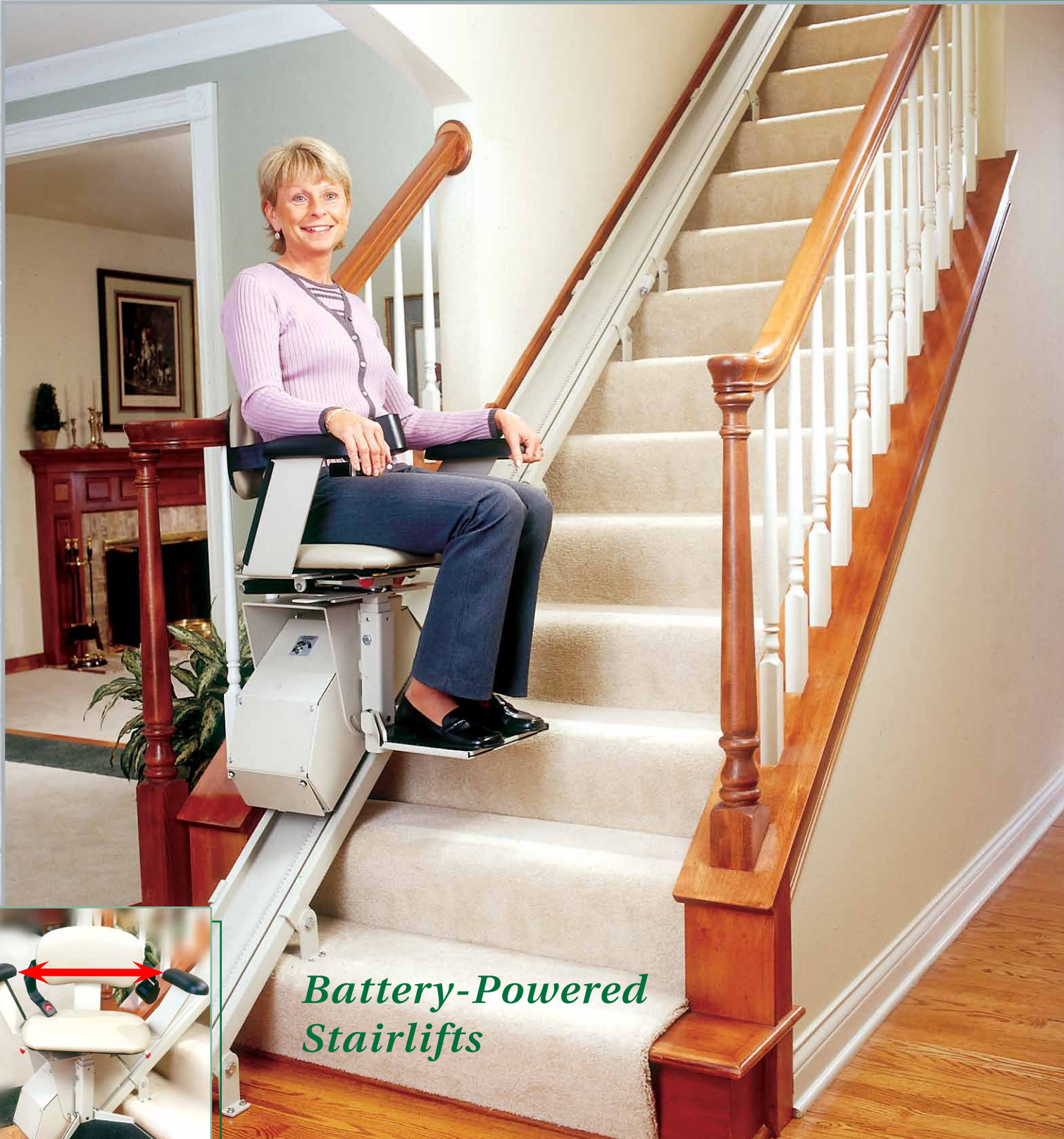

*Battery-Powered  
Stairlifts*



**Adjustable Seat Width**

**BRUNO**  
for your independence

# BRUNO

*The Electra-Ride II Straight Rail Stairlift  
makes going up & down stairs a pleasure!*

## Battery Powered & Adjustable

### **Best Value**

The Electra-Ride II features a wide range of standard features for the best value in stairlifts today! Remote call/sends, 45- and 90-degree seat swivel (at both the top and bottom of the stairs), adjustable swing-out arms and a 350 lb (159 kg) rider capacity are *all* included with the Electra-Ride II. Similarly priced units either don't offer these features or make them available only as extra-cost options.

### **Flexible Seating**

With the Electra-Ride II, seating can be adjusted to meet a wide variety of needs. The seat height can be adjusted from 22.5 in to 25.5 in (572 mm to 648 mm). Plus, the armrests independently swing out for easy wheelchair transfer.

### **Quiet, Powerful, Smooth Ride**

The Electra-Ride II features a powerful rack-and-pinion drive system to give you an ultra-smooth ride. The Electra-Ride II's battery-powered system runs 'whisper quiet' in your home.

### **Reliable Operation**

The Electra-Ride II uses two 12-volt batteries that are powered by a small, unobtrusive battery charger that plugs into any household outlet.

Bruno's leading-edge technology assures you access up and down the stairway *even* if the electrical service in your home is interrupted.



The arms swing out of the way for easy wheelchair transfers using the standard SRE-1550 seat.

The arm release lever is conveniently located under the seat.

Optional, large SRE-1540 seat with fixed arms provides additional seat depth and width.

# BRUNO Electra-Ride™ II

# Compact • Safe Affordable

*With Bruno, your lifestyle doesn't have to be limited by stairs.  
Enjoy the home you love...upstairs and down!*

## Safety

When traveling up and down your stairs, the Electra-Ride II will automatically stop when it "senses" any obstacles. Plus, when you're ready to get off the unit upstairs, the seat easily swivels so you face the hall (a seat belt is standard). The Electra-Ride II meets all ANSI A-18.1 residential requirements as well as all public building specifications (with the optional commercial package).

## Affordable

The Electra-Ride II is simpler and less costly to install than other units since there is no need to install additional electrical outlets...it's completely battery powered!



The arms, seat and footrest flip up, creating plenty of space for guests to walk up the stairs.



The two wireless call/send controls make the installation simple and clean with no wires running along the wall.

*Easy on-and-off thanks to the Electra-Ride II's unique design.*

## Electra-Ride II SRE-1550 Specifications

### Exceptional Value

- ◆ Locking swivel seat, at the top and bottom of the stairway.
- ◆ Two remote, wireless call/send controls.
- ◆ Adjustable armrest width.
- ◆ Footplate safety sensor and seat belt.
- ◆ Charging station indicator lights at the top and bottom of the stairway.
- ◆ On board diagnostics make service easier and more precise.
- ◆ Meets all ANSI - A18.1 requirements.
- ◆ CSA/C-US listed.
- ◆ Bruno's Gold Warranty – Five years for major components and two years on parts.

### Convenient Installation

- ◆ Rail installs to within 10 in (254 mm) of the wall.
- ◆ Can be installed on either side of the stairway with simple modification.
- ◆ Standard track lengths of:
  - 16 ft (4.9 m) – two 8 ft (2.4 m) sections
  - 20 ft (6 m) – two 10 ft (3 m) sections
 (Custom lengths also available)

### Smooth, Powerful Drive System

- ◆ Self-locking worm gear, rack and pinion drive.
- ◆ Maximum rider weight of 350 lbs (159 kgs)
- ◆ Soft start, solid-state controller makes it possible to adjust the speed to the comfort of the user.

### Battery-Powered Technology

- ◆ Ensures access even during power outages.
- ◆ Powered by two 12-volt batteries; continuously charged by a 2-Amp battery charger that plugs into any 110-volt household outlet.

### Comfort

- ◆ Adjustable seat height for easier entry/exit at the top and bottom of the stairway.
- ◆ Swing-away arms for easy wheelchair transfers.
- ◆ Seat is contoured to provide optimum comfort.

BRUNO  
**Electra-Ride™ II**

# BRUNO INDEPENDENT LIVING AIDS, INC.®


*manufactures independence products that help people live their lives to the fullest!*

Ask us about our other...

## INDOOR AND OUTDOOR STRAIGHT RAIL and CURVED RAIL *Stairlifts:*



**ELECTRA-RIDE™ LT** Model SRE-2750

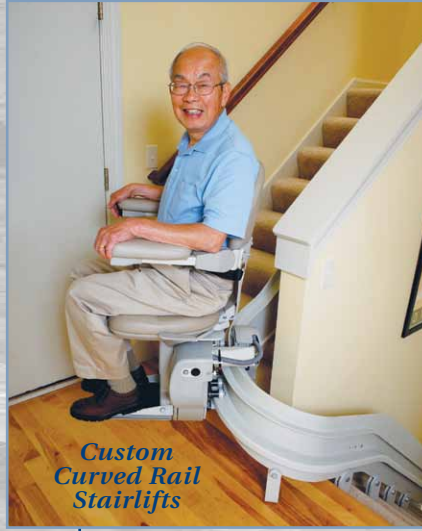

**ELECTRA-RIDE™ III** Model CRE-2110  
The only 400 lb/ 181.8 kg capacity curved stairlift... an **INDUSTRY EXCLUSIVE!**



**OUTDOOR ELECTRA-RIDE™** Model SRE-2010E  
The only 400 lb/ 181.8 kg capacity outdoor stairlift... an **INDUSTRY EXCLUSIVE!**

## OVER 20 DIFFERENT *Vehicle Lifts* INCLUDING:



**OUT-SIDER® MERIDIAN™ SERIES**



**JOEY™ INTERIOR PLATFORM LIFT** VSL-4000



**CURB-SIDER® SERIES**



Accredited by ANSI-RAB NAP  
Accredited by the Dutch Council for Accreditation (DVA)



Founded over 25 years ago, Bruno Independent Living Aids, Inc., is a world-wide manufacturer of accessibility products designed to enhance the lives of those challenged by limited mobility. In pursuing this goal, Bruno is proud to be the only U.S. manufacturer of indoor/outdoor stairlifts, interior/exterior scooter, powerchair and wheelchair vehicle lifts and Turning Automotive Seating™ (TAS™) products to be ISO 9001 Certified. This is one of the most stringent, widely-recognized quality standards in the world and your assurance that Bruno is committed to providing you with the highest-quality, most reliable products in the industry. Ask your certified Bruno dealer for details on all Bruno independence products.

Represented by:

# BRUNO INDEPENDENT LIVING AIDS, INC.®

1780 Executive Drive, P.O. Box 84, Oconomowoc, WI 53066  
www.bruno.com



**TURNY™ ORBIT® SEAT**

All illustrations and specifications in this brochure are based on the latest product information available at the time of publication. Bruno Independent Living Aids, Inc. reserves the right to make changes at any time without notice. Electra-Ride II Brochure 7/10-9 © Bruno Independent Living Aids, Inc. 2010.

Part#: MKT-SRE-1550

